



ESD

Endoscopic Submucosal
Dissection

Minimally invasive endoscopic therapy has been developing and thriving since the Endoscopic Submucosal Dissection (**ESD**) was established more than 10 years ago. The treatment range covers mucous, submucosa, muscularis propria, serosa and even extraluminal lesions.

At the same time, **POEM** (Per Oral Endoscopic Myotomy),



POEM



ESD

STER (Submucosal Tunneling Endoscopic Resection)



STER


And **NOTES** (Natural Orifice Transluminal Endoscopic Surgery) technologies etc. have also been developing rapidly.



NOTES

Endoaccess offers innovative instruments for the minimally invasive endoscopic procedures like ERCP, EMR, ESD, POEM, STER, NOTES etc.

Disposable Ball Type (Ball Type)

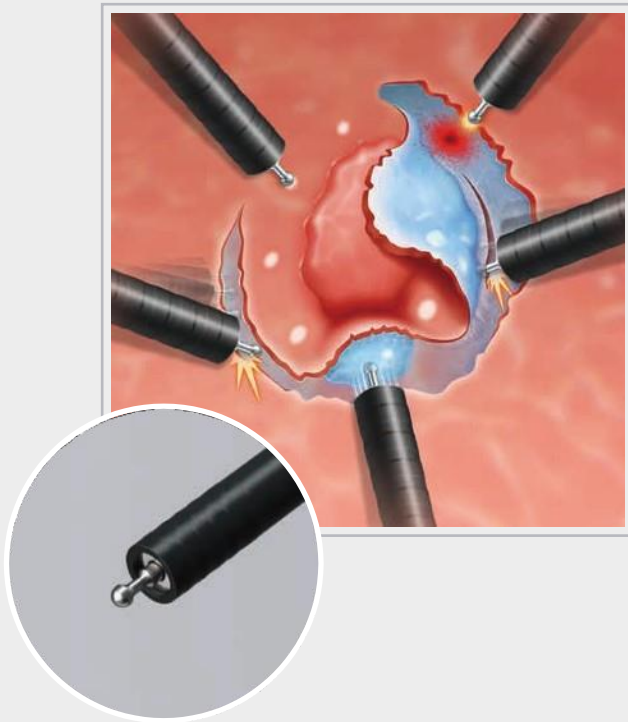

Multifunctional Design:
marking, cutting, dissection, coagulation, spraying, etc.

Ball Shape Type:
Good operability, precise coagulation effect. Good traction during cutting and peeling

Spray function:
For spraying into the submucosa to lift the mucosa or to show the bleeding area

Product Specifications

Item no	Tube \varnothing and working length	Working channel \varnothing	Application
EA-S24BT18	\varnothing 2.4 x 1800mm	\geq 2.8mm	Gastroscopy
EA-S24BT23	\varnothing 2.4 x 2300mm	\geq 2.8mm	Colonoscopy

Instructions



Marking Coagulation (tip in) Cutting Dissection (tip out) Spray (tip out)

Disposable Ball Type (CT Knife)



Insulated Tip Design:

Ensures safety when cutting and dissecting
 Excellent conductivity
 Gentle cutting and stable coagulation

Adjustable electrode length:

Adjustable from 0 mm to 4 mm according to operation procedure and lesions


Spray function:

For spraying into the submucosa to lift the mucosa or to show the bleeding area


Product Specifications

Item no	Tube Ø and working length	Working channel Ø	Application
EA-S24B3T18	Ø 2.4 x 1800mm	≥ 2.8mm	Gastroscopy
EA-S24B3T23	Ø 2.4 x 2300mm	≥ 2.8mm	Colonoscopy

Instructions

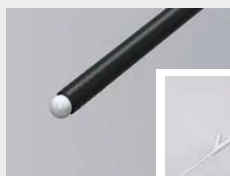

Tip in

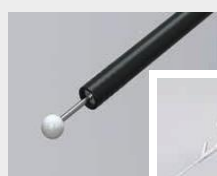

Tip out



Spray (tip out)



Disposable Cylindrical Type (Tunnel Knife)



Multifunctional Design:

Marking, cutting, injection, lifting, dissection, coagulation, etc.

Needle shape tip:

Sharp cutting tip extendable up to 2 mm

Submucosa injection function:

High quality hydraulic solution pressure with needle point design; the knife can inject precisely into the submucosa and achieve a perfect mucosal lifting effect

Product Specifications

Item no	Tube \varnothing and working length	Working channel \varnothing	Application
EA-S24ST18	\varnothing 2.4 x 1800mm	\geq 2.8mm	Gastroscopy
EA-S24ST23	\varnothing 2.4 x 2300mm	\geq 2.8mm	Colonoscopy

Instructions




Marking /coagulation (tip in)

Cutting /dissection (tip out)

Spray (tip out)

Disposable Cylindrical Type (Square Knife)

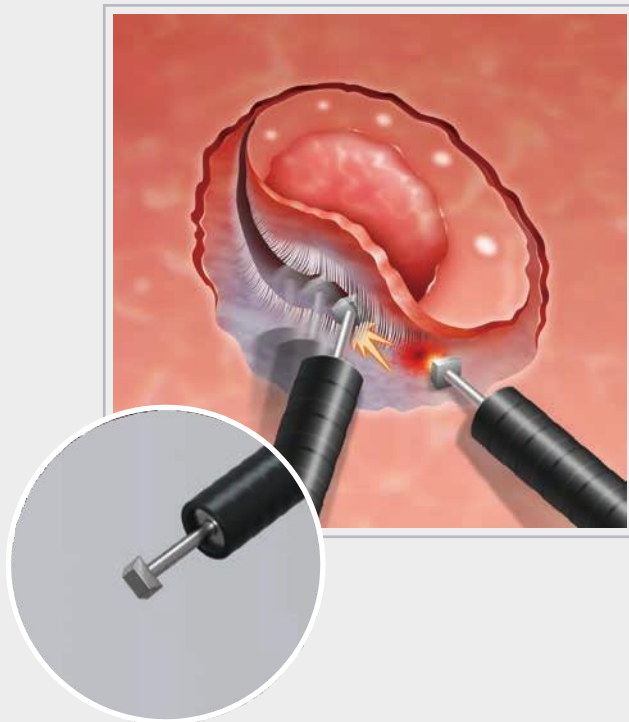

Square Shape Tip Design:

Marking and coagulation
The design of the tip allows to lift or hook the tissue to cut efficiently. Good support surface and better coagulation effect

Adjustable electrode length:

Adjustable cutting tip from 0 to 4mm


Spray function:

For spraying into the submucosa to lift the mucosa or to show the bleeding area effect

Product Specifications

Item no	Tube Ø and working length	Working channel Ø	Application
EA-S24S2T18	Ø 2.4 x 1800mm	≥ 2.8mm	Gastroscopy
EA-S24S2T23	Ø 2.4 x 2300mm	≥ 2.8mm	Colonoscopy


Instructions

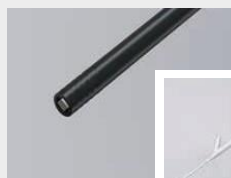

Marking /coagulation (tip in)



Cutting /dissection (tip out)

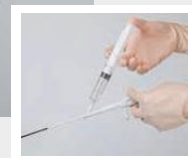

Inject (tip out)



Endoaccess GmbH
Feldriethe 1
D – 30826 Garbsen
Germany
info@endoaccess.de
www.endoaccess.de